

# BRODA<sup>®</sup>

# Comfort

---

## Rehab Wheelchair



HCPCS  
E1161

*Shown with optional head support, Matrx Elite back support, Matrx Vi seat cushion, and angle-adjustable lower leg supports*

## **Patient, Facility, & Caregiver Benefits**

---

- Pressure redistribution that aids in maintaining skin integrity
- Durable frame to withstand heavy workload healthcare environments
- Adaptable for wide range of patients and care needs



*Front pivot seat tilt maintains user comfort with smooth, effortless tilt*



*Posterior tilt optimizes pressure redistribution*



*Solid seat pan offers support and is compatible with a variety of seat cushion options*



*Range of simple adjustments accommodates a variety of needs*

**brodaseating.com | 1-844-552-7632**

Matrix® seating options available to ship preconfigured on-chair.

## matrix® SEATING SERIES

### Loxx Headrest



Highly adjustable hardware with tool-less adjustment levers for easy fitting and position changes.

### E2 Back



**Shell** - Contoured aluminum with cutouts for lifting and breathability. Available in standard (3") OR deep (6") contours.  
**Foam** - HR (high resilient) foam enclosed in cover.

### Libra Cushion



**Foam** - HR (high resilient) molded foam made with UltraFresh antimicrobial for odor protection. Contoured pelvic well, waterfall front edge.



Additional accessories available: lower leg supports, trays, and positioning belts

## SIMPLIFIED DURABILITY FOR REHABILITATION CARE

### SPECIFICATIONS

<b>Degrees of Tilt</b>	Up to 45°
<b>Weight Capacity</b>	Up to 250 lbs. (114 kg)
<b>Angle Adjustable Back Canes</b>	5°, 10°, 15°, 20°
<b>Seat Widths</b>	16", 18", 20", 22" (40cm, 45cm, 50cm, 55cm)
<b>Seat Depths</b>	Adjustable: 16", 17", 18", 19", 20" (40cm, 43cm, 45cm, 48cm, 50cm)
<b>Seat to Surface Heights</b>	Adjustable: 14", 15", 15.5", 16", 16.5", 17", 17.5", 18", 19" (35cm, 38cm, 39cm, 42cm, 43cm, 44cm, 45cm, 48cm)
<b>Wheels</b>	5", 6" or 8" Front Caster Wheels with fork (13cm, 15cm, 20cm) 20", 22" or 24" Rear Mag Wheels (50cm, 55cm, 60cm)
<b>Lifetime Warranty Frame</b>	10 years on 16 Gauge tubular steel frame
<b>Government Coding (US &amp; CA)</b>	HCPCS: E1161 ADP (ON): WMA5B-025

### FRAME OPTIONS

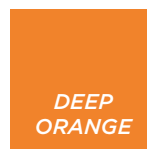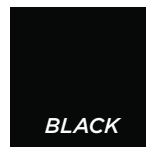

**brodaseating.com**  
**1-844-552-7632**