

Synthesis

Rehab Wheelchair

REHAB
BRODA[®]



Tilt functionality provides additional pressure offloading



Front pivot seat tilt maintains a forward line of sight, enhancing engagement.



Adjustable upper lateral supports fit a variety of users



Solid seat pan provides support for a variety of cushions

Patient, Facility, & Caregiver Benefits

- Increases quality of life and engagement with surroundings
- Promotes safe patient handling with enhanced positioning
- Proper seating alignment reduces patient fatigue and discomfort

matrx[®] SEATING SERIES

The Synthesis of Comfort, Functionality, and Style.

SPECIFICATIONS



Libra Cushion

Foam - HR (high resilient) molded foam made with UltraFresh antimicrobial for odor protection. Contoured pelvic well, waterfall front edge.

Vi Cushion



Ischial Insert - Supersoft viscoelastic foam insert for envelopment and protection of bony prominences.



PSP Cushion



Polymer Insert - Molded-in polymer gel provides extra skin protection and shear reduction

<i>Infinitely Adjustable Tilt</i>	<i>up to 40°</i>
<i>Infinitely Adjustable Recline</i>	<i>up to 90°</i>
<i>Weight Capacity</i>	<i>350 lbs. (158 kg)</i>
<i>Seat Widths</i>	<i>16", 18", 20", 22", 24", 26" (40cm, 45cm, 50cm, 55cm, 61cm, 66cm)</i>
<i>Seat Depths</i>	<i>17", 18.5", 20", 22" (43cm, 46cm, 50cm, 56cm)</i>
<i>Seat Surface Heights</i>	<i>19.5", 20.5", 22", 23.5" (48cm, 51cm, 55cm, 59cm)</i>
<i>Wheels - Mag Wheels</i>	<i>14", 20" Rear mount mag wheels (35cm, 50cm, 55cm) 14" Mid mount mag wheels (35cm)</i>
<i>Wheels - Caster Wheels</i>	<i>5" Front total lock or non-locking caster wheels (13cm) 6" Rear multi-directional locking caster wheels (15cm)</i>
<i>Lifetime Warranty Frame</i>	<i>10 years on 16 Gauge tubular steel frame</i>
<i>Pressure-Reducing Seat & Back Support</i>	<i>Comfort Tension Seating[®] Back & Seat Supports Optional Matrx[®] Seating Series Cushions</i>
<i>Standard Padding</i>	<i>Neck support, upper lateral supports, lower lateral supports, and sole, foot, seat, and back cushions</i>

FRAME OPTIONS

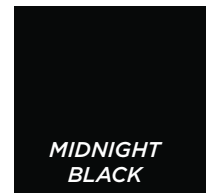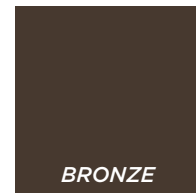

Additional accessories available: foot supports, specialty padding packages, trays, and positioning straps